

Jannali

Multi-function Bed

ALIVIO

Enhanced comfort, convenience, safety, and support



Soft fabric (Black) shown



The **Jannali multi-function bed** offers a range of benefits for those seeking enhanced comfort, convenience, safety, and support in their sleep or rest routines.

Jannali's customisable positioning features enable you to find your optimal sleeping or resting position – with five (5) different pre-set positions and two memory functions using the bluetooth phone App (iOS & Android compatible) to store your preferred bed positions.

An Autoregression section is built into the bed base to decrease pressure on the abdominal area and reduce shear forces on the occupant. The Jannali also features Trendelenburg, and reverse Trendelenburg positioning.

The electric height adjustment function provides independence and reduces reliance on carers for assistance. For occupants requiring assisted transfers it provides the carer with flexibility of height depending on the transfer equipment being used. For additional safety, a non-entrapment zone under the bed stops the bed frame from touching the floor. The standard under bed dimming LED light (motion detection version is optional) provides visibility at night.

The full-featured Jannali has various convenient features, including a standard wired USB handset (or downloadable phone App), battery backup, and USB ports for charging devices.

Optional: The Jannali also comes with an exciting range of options, including matching headboard/footboard, side/corner mattress retainers, armrest/armrest cover, short/long bedrails, side pocket, height adjustable bed stick, height adjustable assist rail and motion detection LED light.

Jannali comes in a range of modern, stylish colours and fabrics, including a soft-touch easy to clean brushed vinyl (Black) and soft leather-look PU option (Beige and Canley), with colour availability for optional matching companion beds.

Fabric/Colour Options

Black

Canley

Beige

Spare Parts: a range of spare parts are available. Please consult with your Jannali dealer.

Warranty:

Bed frame _____ 10 years

Actuator and control box _____ 5 years

Wireless remote-control handset _____ 1 year


The Jannali bed will suit most standard and pressure-care mattresses – please consult your Jannali dealer before purchasing.

Functional Features

- Electrically adjustable height, Trendelenburg and reverse Trendelenburg, backrest and leg rest profiling
- Auto Regression backrest reduces abdominal pressure and shear forces
- A non-entrapment zone stops the bed from touching the floor for additional safety
- Wired USB handset
- Multiple One Touch Buttons for simple operation (Zero Gravity, TV, Relax, Flat, Home)
- Memory setting to user's favourite position
- Home function to bring the bed to lowest and flat position
- Under bed dimming LED light as standard (motion detection version is optional)
- 6 USB ports – 4 x USB-A , 2 x USB-C
- Backup battery
- Control from your smart phone via Bluetooth and free App (Apple and Android)
- Swivel and lockable castors

Options


- Short/Long bed rails (1)
- Height adjustable assist rail. Easy bed positioning and removable (2)
- Height adjustable bed stick, positionable and removable
- Headboard/footboard (3)
- Wireless remote-control handset
- Simplified grey wireless remote-control handset (KA3200RCS) (4)
- Side/corner mattress retainers
- Motion detection LED light
- Matching companion bed
- Three Color/Fabric Options – soft fabric (Black) and soft leather-look PU (Beige and Canley)

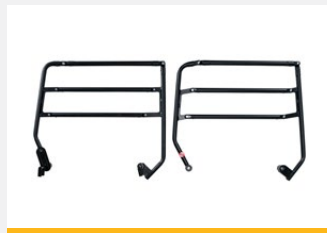

1. Long bed rails



2. Height adjustable assist rail

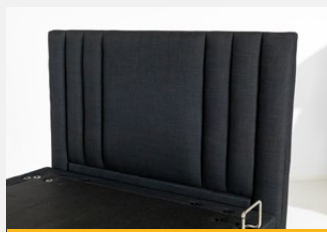

3. Headboard



4. Simplified wireless remote-control handset

	Long Single	King Single
Product weight without Accessories	100kg	108kg
Width	920mm	1060mm
Height (floor to top of mattress support platform)	220mm to 750mm	
Length	2030mm	
Safe Working Load (SWL)	200kg (User weight+ mattress weight+ accessories weight)	
Max User Weight (MUW)	= SWL- Mattress - all accessories installed. 160kg Typical	
Trendelenburg and Reverse Trendelenburg	0-14 degree	
Backrest angle	0-75 degree +/- 3 degree	
Leg rest angle	0-45 degree +/- 3 degree	
Castors	8 swivel locking castors 75mm / 3"	
Input Voltage	100-240V	
Output Voltage	DC29V	
Current	3.5A	
Power consumption	50W	
Code	Black - KA3200LSBLK, Canley - KA3200LSCAN Beige - KA3200LSBEI	Black - KA3200KSBLK, Canley - KA3200KSCAN, Beige - KA3200KSBEI
Companion Code	Black - KA3200LSCBLK, Beige - KA3200LSCBEI Canley - KA3200LSCCAN	Black - KA3200KSCBLK, Beige - KA3200KSCBEI Canley - KA3200KSCCAN

Please contact customerservice@kcare.com.au for more information.

www.kcare.com.au | 1300 783 783